

Analyst

Accepted Manuscript

This article can be cited before page numbers have been issued, to do this please use: S. Garcia, J. Herniman, P. Birkin, J. Pilgrim, P. Lewin, G. Wilson, J. Langley and R. C. D. Brown, *Analyst*, 2020, DOI: 10.1039/D0AN00602E.



This is an Accepted Manuscript, which has been through the Royal Society of Chemistry peer review process and has been accepted for publication.

Accepted Manuscripts are published online shortly after acceptance, before technical editing, formatting and proof reading. Using this free service, authors can make their results available to the community, in citable form, before we publish the edited article. We will replace this Accepted Manuscript with the edited and formatted Advance Article as soon as it is available.

You can find more information about Accepted Manuscripts in the [Information for Authors](#).

Please note that technical editing may introduce minor changes to the text and/or graphics, which may alter content. The journal's standard [Terms & Conditions](#) and the [Ethical guidelines](#) still apply. In no event shall the Royal Society of Chemistry be held responsible for any errors or omissions in this Accepted Manuscript or any consequences arising from the use of any information it contains.

ARTICLE

Quantitative UHPSFC-MS analysis of elemental sulfur in mineral oil via derivatization with triphenylphosphine: Application to corrosive sulfur-related power transformer failureReceived 00th January 20xx,
Accepted 00th January 20xx

DOI: 10.1039/x0xx00000x

Sérgio B. Garcia,^a Julie Herniman,^a Peter Birkin,^a James Pilgrim,^b Paul Lewin,^b Gordon Wilson,^c G. John Langley^{*a} and Richard C. D. Brown^{*a}

An ultrahigh-performance supercritical fluid chromatography–mass spectrometry (UHPSFC-MS) method has been developed as a rapid and reliable analytical method for the detection and quantification of elemental sulfur in mineral transformer oil. The method described in this paper is based on the selective reaction of elemental sulfur with triphenylphosphine (TPP). The derivatization of the elemental sulfur requires minimal sample preparation and resulted in the formation of a single compound, namely triphenylphosphine sulfide (TPPS). This derivative is quantified from the complex oil composition using electrospray ionization-mass spectrometry (ESI-MS) in selected ion monitoring (SIM) mode, and the reported UHPSFC-MS method allows detection and quantification of the derivative at ppb levels. As sulfur contamination in mineral transformer oil has been linked to costly failures of oil/paper-based power transformers due to corrosion, the analytical approach is demonstrated through its application to mineral oil samples from in service and decommissioned power transformers. The method is ideal as a routine test or to confirm the presence of elemental sulfur in samples where corrosion has occurred.

Introduction

Oil insulation systems are used by numerous power apparatus. Key components in power networks, like power transformers, use mineral oil not only for cooling but also for electrical insulation. High voltage (HV), super grid transformers cost millions of pounds. Therefore, reliability and economic viability of power distribution is highly dependent on availability, performance and longevity of power transformers.¹

In the last decades, the risk of failure in electrical equipment owing to the presence of corrosive sulfur in the mineral oil has increased. This risk has been understood for many decades but actual failures world-wide were very rare until around 2000. The International Council on Large Electric Systems (CIGRÉ) estimated the total number of failures due to corrosive sulfur to be in the region of 100 between 2000 and 2009.² A number of the cases were associated with the presence of silver sulfide in silver-coated contacts, mainly in the tap-changers.³ Formation of surface metal sulfide will lead to overheating of the contacts followed by the detachment of conductive particles. These particles can cause breakdown and rated voltage stress.²

Transformer failures and/or repairs are expensive and

should be avoided whenever possible. Outages if not predicted in advance can be expensive, not only for generators but also for networks. To avoid unexpected failures, diagnostic measures are necessary. Standard methods for corrosive sulfur detection in mineral oil have been published.^{4–8} These methods are not specific for S₈, they are time consuming and subjective due to a visual rating process. Hence, analytical methods that allow rapid, reliable and reproducible quantification of S₈ in mineral transformer oil are needed.

Several analytical techniques have been applied to determine S₈ in hydrocarbons, but none specifically in mineral insulating oil. The use of gas chromatography (GC) and a number of variants (e.g., use of a selective detector for sulfur) can be found in the literature.^{9–13} Other methods include the use of polarography,^{14–15} liquid chromatography (LC),¹⁶ and GC-Inductively Coupled Plasma - MS (GC-ICP-MS)^{17,18} There are several issues that can arise when applying these methods to quantify S₈ in transformer mineral oil.

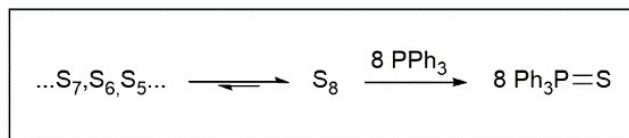


Fig. 1 Derivatization reaction.

In mineral oil, at room temperature elemental sulfur exists in an allotropic equilibrium ($\dots S_6 \rightleftharpoons S_7 \rightleftharpoons S_8 \dots$) with S₈ being the predominant allotrope,¹⁹ at elevated temperatures, as present in a standard GC injector, the allotropic equilibrium shifts to

^a School of Chemistry, University of Southampton, Southampton SO17 1BJ, United Kingdom. *E-mail: rcb1@soton.ac.uk; G.J.Langley@soton.ac.uk; Tel: +44 (0)23 80594108.

^b The Tony Davies High Voltage Laboratory, University of Southampton, Southampton SO17 1BJ, United Kingdom.

^c National Grid, National Grid House, Warwick Technology Park, Gallows Hill, Warwick CV34 6DA, United Kingdom.

lower mass allotropes and this ensures that elemental sulfur does not elute from GC columns as a single and narrow chromatographic peak. Mineral oil has miscibility issues with reversed phase HPLC solvents, while normal phase HPLC is more challenging to interface to MS. Even coupled with GC, ICP-MS systems would suffer from matrix interference due to the myriad hydrocarbons present in mineral oil.²⁰ UHPSFC can overcome the solubility issues related to the mineral oil, and when combined with a suitable derivatisation method and ESI-MS detection, may deliver simplified detection with minimal interference from the matrix.

One derivatisation method, consists of the reaction of sulfur, independent of its allotropic composition, to give a single stable compound, namely triphenylphosphine sulfide (TPPS, Fig. 1).²¹ This derivatisation method has been employed to determine elemental sulfur in different media.^{16,22} Another clear advantage of the derivatisation is an almost 9-fold increase in the mass concentration of the sulfur containing species (1 µg/mL of sulfur → 9.18 µg/mL of TPPS).²¹

Herein, we describe a novel UHPSFC-MS method via derivatisation with TPP that allows the detection of elemental sulfur in mineral oil down to ppb levels. The reported method is capable of maintaining a high sensitivity and accuracy when applied to real mineral transformer samples.

Materials

Chemicals

Toluene (High-performance liquid chromatography (HPLC) grade) and methanol (Liquid Chromatography – Mass Spectrometry (LC-MS) grade) were purchased from Fisher Scientific (Loughborough, UK). Food grade carbon dioxide was purchased from BOC Special Gases. Elemental sulfur (99.998%), TPP (99 %), TPPS (98 %) and ammonium acetate were purchased from Sigma–Aldrich (Gillingham, U.K.). Formic acid was purchased from Riedel-de Haën (Seelze, Germany) and National Grid plc (London, UK) kindly supplied representative mineral oil samples.

Mineral oil samples

National Grid plc provided mineral oil samples for analysis, some of which are believed to possess corrosive properties and all with varying compositions and chemistries. The samples were collected from a diverse range of transformers, presenting different working conditions and ages.

Experimental

Analytical procedure

Stock solution and calibration

Table 1 Inlet method - UHPSFC

Time (min)	Solvent A (CO ₂) %	Solvent B (methanol + 25 mM ammonium acetate) %
0.00	98	2
2.00	60	40
2.30	60	40
3.00	98	2

A working stock solution of TPPS was made in 50:50 mineral oil/toluene (v/v) at a concentration of 100 µg/mL. Dilutions of the working solution with 50:50 mineral oil/toluene (v/v) were used to prepare calibration standards. The calibration standards were prepared by dilution with 50:50 mineral oil/toluene (v/v) to give concentrations of 0, 0.05, 0.10, 0.25, 0.50, 1.0, 2.5 and 5.0 µg/mL of TPPS. These provide concentrations at the lower limit of the quantification (LLOQ) and upper limit of quantification (ULOQ) for TPPS.

Sample preparation and derivatization procedure


Mineral oil samples were stirred at room temperature, and then 5.0 mL of each sample was transferred to a round bottom flask to which 5.0 mL of toluene was added. For derivatization, 20 mg of TPP was added to each round bottom flask and the mixture was stirred at 60 °C. Sufficient TPP is available for oil samples with a concentration of elemental sulfur up to 490 µg/mL. In this work, the samples were allowed to react for 3 h before being diluted 1:100 with 50:50 new mineral oil/toluene(v/v) and loaded into autosampler vials. All samples were analysed by UHPSFC-MS immediately after preparation.

UHPSFC method

An Acquity UPC² system (Waters Corp., Milford, MA, USA) was employed for all analyses. This system is equipped with a binary solvent manager, a heated column manager and an internal autosampler with a cooled sample tray. A HSS C18 SB column (1.8 µm, 3.0 x 100 mm) was used. The column was held at 40 °C in a column oven and 2.0 µL of each sample solution was injected. A mobile phase modifier, methanol (25 mM ammonium acetate) was used. The supercritical (sc) CO₂ back pressure of the system was set to 150 bar. A 3 min gradient elution was performed using the method described in Table 1 and a 1 min pre-run was used for column equilibration.

Before the analyses, three blank injections of methanol were performed to ensure system stability and cleanliness. Analyses commenced with triplicate injections of the calibration curve solutions (low concentration to high), each standard solution was followed by two blank injections. Study samples (triplicate measurements) were then analysed in a randomised order to minimise bias due to batch effects. The sequence of analysis is shown in Fig 2.



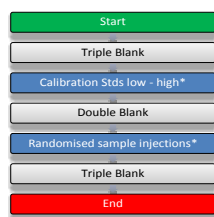

* A blank was conducted in between samples.

Fig 2 Sequence of the samples during the analyses

ESI-MS method

A Waters TQD (Waters Corp., Milford, MA, USA) triple quadrupole mass spectrometer was used. An Internal Solvent Manager (ISM) was used to introduce make up solvent, methanol + 1% formic acid, after chromatographic separation at a flow rate of 0.5 mL/min to enhance ionisation and ensure a stable spray. The operating conditions for the positive ion ESI were as follows: capillary voltage 3.5 kV, cone voltage 20 V, source temperature 150 °C, desolvation temperature 500 °C. Nitrogen was used as the desolvation gas at a flow of 300 L/h. Full scan data over the m/z range of 100 – 800 were acquired for each injection. A positive ion ESI mass spectrometry single ion monitoring (SIM) method was created for the ion of interest m/z 295, which corresponds to protonated TPPS [TPPS + H]⁺.

Data processing

Waters MassLynx™ software (version 4.1) (Waters Corp., Milford, MA, USA) was used to collect and process data for all the UHPSFC-MS experiments.

Results and discussion

Method development

TPP reacts quantitatively with elemental sulfur to give TPPS,²⁷ which displays favourable ionization properties under positive ESI conditions. Potential considerations associated with the use of TPP as derivatisation reagent are the possible occurrence of side reactions, such as oxidation to triphenylphosphine oxide (TPPO). However, probably due to the presence of the anti-oxidant (2,6-di-tert-butyl-*p*-cresol) in mineral oils, oxidation was not observed, and no modification of the derivatisation reaction conditions was required.

Different chromatography-MS approaches were investigated for quantification of TPPS in mineral oil. Electron ionization (EI) provided an efficient method for detection of TPPS, however GC-MS proved to be unsuitable for reliable analysis of samples when mineral oil was present due to ion source contamination and subsequent rapid loss of sensitivity.

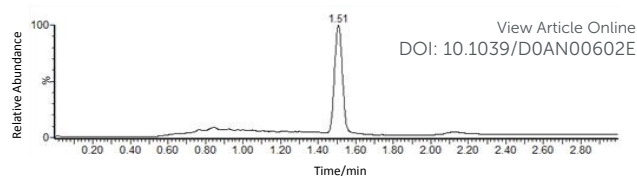


Fig 3 ESI UHPSFC-MS SIM of TPPS at 1 µg/mL (S/N = 161).

Therefore, atmospheric pressure ionisation techniques were explored in combination with SFC to provide reliable and short analysis times. The effectiveness of the ionization of TPPS was explored using ESI for a 1 µg/mL standard solution, using a Single Ion Monitoring (SIM) method (m/z 295). TPPS was found to undergo efficient ionisation, with reduced interference arising due to ionisation of the matrix. Under the selected conditions TPPS has a retention time of 1.5 min, and the protonated molecule ([TPPS + H]⁺) is observed with a S/N of 161 (Fig 3).

Positive ESI MS conditions for TPPS and the chromatographic method were optimized providing an analysis time of 3 min/sample. For quantification, an external standard calibration was employed where the ion at m/z 295 ([TPPS + H]⁺) was used. Sample preparation was adapted to mitigate against the low solubility of TPPS in neat mineral oil, different solvents were studied based on their capacity to dissolve both TPPS and mineral oil. A 50:50 mineral oil/toluene mixture (v/v) provided acceptable solubility while maintaining good chromatographic efficiency and resolution.

Method validation

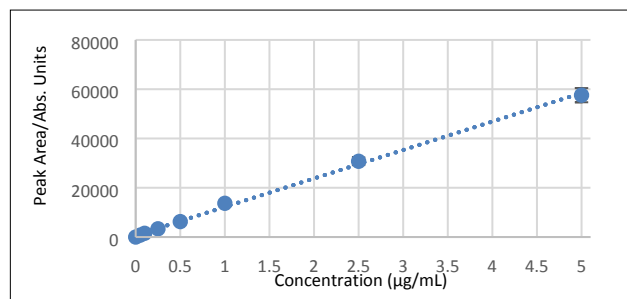


Fig. 4 Calibration curve for TPPS (number of replicate measurements = 4).

A range of analytical attributes were investigated, including linearity, lower and upper limits of quantification (LLOQ and ULOQ), reproducibility, carryover, matrix interference and derivative stability. Standard curves were found to be linear over the range measured 0 – 5 µg/mL, with correlation coefficients (R^2) of 0.992 or better as determined over the three days used for validation. A representative calibration curve is shown in Fig 4. A LLOQ of 0.05 µg/mL (S/N of 33) and an ULOQ of 5 µg/mL was used for all analyses. The concentration range covered by the standard curves equates to 0 – 10 µg/mL of TPPS due to the 50:50 dilution of the samples prior to the derivatisation reaction. This concentration of TPPS corresponds to a range of 0 – 1.09 µg/mL of elemental sulfur. Carryover was assessed by calculating the concentration of TPPS present in the



blanks performed in between standard analysis, and was only found present when analysing samples with concentrations close to the ULOQ and presented a value of $\leq 0.5\%$. TPPS was found to be stable for at least 24 h at ambient temperature or in the autosampler at 4 °C. Intraday reproducibility was determined by the analysis of three standard samples with different concentrations over the range of 0 - 5 $\mu\text{g}/\text{mL}$ on three separate occasions, the method was found to be reproducible with a maximum relative standard deviation (RSD) of 2.03%. In this work, a new calibration curve was prepared daily. Interference of the matrix on the concentration of TPPS is minimised as most components in mineral oil are not ionized by ESI, and a SIM method is employed. Furthermore, the retention time of TPPS is separated from the vast bulk of the oil components.

Derivatisation selectivity and efficiency

The selective reaction of TPP with elemental sulfur to form TPPS is the basis of this analytical method. Side reactions of the TPP with the other sulfur compounds present in the mineral oil would present false positives/increased concentrations, and are clearly undesirable. TPP reactivity with other representative sulfur compounds was studied using the conditions of the derivatisation method. It is not practical to test all sulfur compounds that are believed/known to be present in the mineral oil, so three representative sulfur compounds with different reactivities were studied. The compounds were dibenzyl disulfide (DBDS), benzothiophene (BT) and 4,6-dimethyldibenzothiophene (4,6-DMDBT). Each was reacted with TPP and analysed using the UHPSFC-MS SIM method. From the results obtained, it can be concluded that there are no products of the reaction between BT or 4,6-DMDBT with TPP. The reaction of the DBDS with TPP, presents only one trace product (TPPS) with a yield of $< 1\%$. This infers that the side reactions of the TPP with typical organosulfur compounds present in the oil are also likely to be negligible, giving confidence that TPPS formed during our studies will arise from the presence of elemental sulfur. These findings were later confirmed by GC-MS and ^{31}P NMR experiments. The progression/conversion of the reaction was studied and determined by the reaction of a known concentration of elemental sulfur with TPP (excess, 1.1 eq.) using the conditions described in the derivatisation procedure. The concentration used was chosen in order to allow the reaction to be monitored by ^{31}P NMR and by our UHPSFC-MS SIM method. After 3 h, integration of the signal of TPPS in the ^{31}P NMR (-4.58 ppm, CDCl_3 , 162 MHz) showed approximately 98% conversion of S_8 to TPPS. The reaction was monitored for two more hours and this value did not change further. At the same time and during the same period, reaction progress was followed using our UHPSFC-MS SIM method, triplicate measurements were conducted after 1, 1.5, 2, 3.5 and 5 h. The results are shown in Fig 5, taking into consideration the results obtained from both experiments the derivatisation reaction time was set to 3 h.

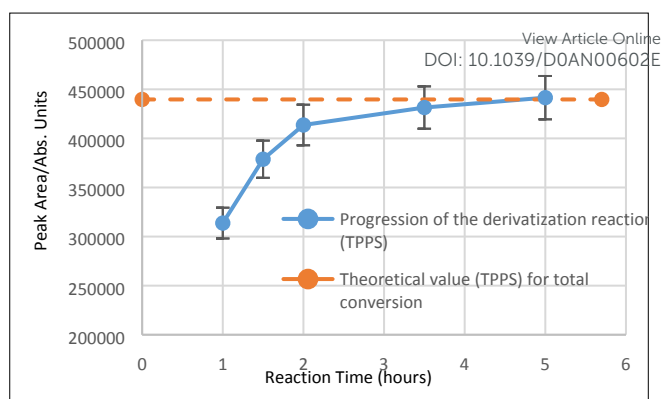


Fig. 5 Efficiency of the derivatisation reaction using the UHPSFC-MS.

Analysis of real mineral oil samples

Following validation, the method was then applied to the analysis of real mineral oil samples from power transformers. To test the reproducibility of this method when applied to real mineral oil samples, two different samples supplied by National Grid were analysed 3 times and on different days (Table 2). The example shown in Fig 6, is the SIM chromatogram obtained from a real mineral sample oil that was doped with elemental sulfur prior to the derivatization process. Inspection of the chromatogram shows a peak at the retention time of 1.62 min,²⁰ corresponding to TPPS, and showing the protonated molecule $[\text{TPPS} + \text{H}]^+$ in the mass spectrum. Two peaks at 2.06 and 2.21 min, elute at retention times separated from TPPS and do not compromise accurate quantification. Furthermore, the different compositions inherent to each sample, due to different working conditions, did not have any negative effect on the analyses performed in this study.

Table 2 Reproducibility of the method when analysing real mineral oil samples

*Values are in $\mu\text{g}/\text{mL}$ of S_8		Sample A	Sample B
Day 1	Experiment 1	0.539	0.739
	Experiment 2	0.542	0.737
Day 2	Experiment 3	0.524	0.694

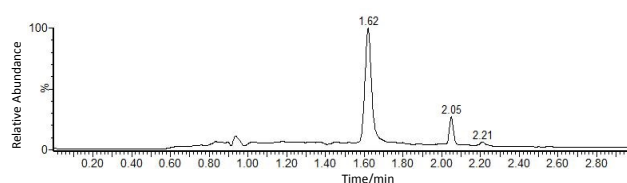


Fig. 6 SIM chromatogram obtained from a real mineral oil sample doped with elemental sulfur (TPPS peak at 1.62 min).



Conclusions

The use of triphenylphosphine as a derivatisation reagent enabled high sensitivity and targeted analysis of elemental sulfur in mineral oil using an UHPSFC-MS method. The derivative, TPPS, enables sensitive, reproducible and specific analysis down to 0.05 µg/mL of TPPS. The use of ESI MS generates the protonated molecule, and results in low interference and background noise from the matrix. This provides a method capable to detect elemental sulfur in samples with varying composition and chemistry, including used/in service oils. This new UHPSFC-MS method has proven to be a fast, robust and reliable tool to identify oil samples that might present a risk for the power transformer.

24 The shift of the retention time of the TPPS is due to small changes in the composition of solvent B.

View Article Online
DOI: 10.1039/D0AN00602E

Conflicts of interest

There are no conflicts to declare.


Acknowledgements

The authors acknowledge financial support from National Grid and EPSRC (EP/K039466/1).

Notes and References

- 1 M. Heathcote, J & P Transformer Book, 2007, **13**, 40.
- 2 CIGRE Technical Brochure n° 378, WG A2.32, 2009, **1**, 5.
- 3 CIGRE Technical Brochure n° 625, WG A2.40, 2016, **1**, 74.
- 4 IEC 60296, Fluids for Electrotechnical Applications: Unused Mineral Insulating Oils in Electrical Equipment, 2012, **4**, 14.
- 5 IEC 60422, Supervision and Maintenance Guide for Mineral Insulating Oils in Electrical Equipment, 2013, **4**, 11.
- 6 Standard DIN 51353, Detection of Corrosive Sulphur, 1985, **1**, 1.
- 7 ASTM D1275-15, Standard Test Method for Corrosive Sulfur in Electrical Insulating Liquids, 2015, **4**, 1.
- 8 ISO 5662, Petroleum products – Electrical Insulating Oils – Detection of Corrosive Sulfur, 1997, **2**, 1.
- 9 L.G. Borchardt and D.B. Easty, *J. Chromatogr.*, 1984, **299**, 471.
- 10 D.L. Struble, *J. Chromatogr. Sci.*, 1972, **10**, 57.
- 11 J.J. Richard, R.D. Vick and G.A. Junk, *Environ. Sci. Technol.*, 1977, **11**, 1084.
- 12 J.L. Guinon, J. Monzo, J. Garcia-Anton, C. Urena and J. Costa, *Fresenius' J. Anal. Chem.*, 1990, **337**, 372.
- 13 P.D. Clark and K.L. Lesage, *J. Chromatogr. Sci.*, 1989, **27**, 259.
- 14 H.S. Katal, A. A. Miran Beigi, M. Farazmand and S. A. Tash, *Analyst*, 2000, **125**, 903.
- 15 B.R. Olofsson, *Anal. Chim. Acta*, 1985, **177**, 167.
- 16 B.F. Taylor, T.A. Hood and L.A. Pope, *J. Microbiological Methods*, 1989, **9**, 221.
- 17 G. D. Woods and F. I. Fryer, *Anal. Bioanal. Chem.*, 2007, **289**, 753.
- 18 W. M. Geiger, S. McSeehey and M. J. Nash, *J. Chrom. Sci.*, 2007, **45**, 677.
- 19 B. Meyer, *Chem. Rev.*, 1976, **76**, 367.
- 20 J. Heilmann and K. G. Heumann, *Anal. Chem.*, 2008, **80**, 1952.
- 21 P.D. Bartlett and G. Meguerian, *J. Amer. Chem. Soc.*, 1956, **78**, 3710.
- 22 R. E. Pauls, *J. Chrom. Sci.*, 2010, **48**, 283.
- 23 ASTM D7800/D7800M-14, Standard Test Method for Determination of elemental sulfur in Natural Gas, 2014, **5**, 1.



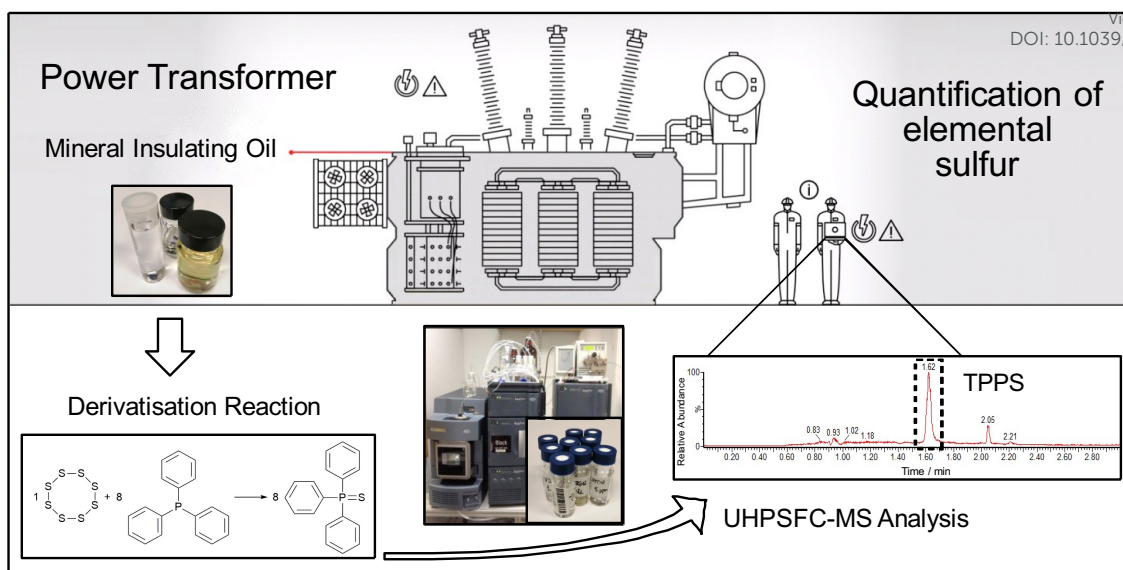

Quantification of elemental sulfur in power transformer mineral oil is achieved using UHPSFC-MS analysis of an easily accessible derivative $\text{Ph}_3\text{P}=\text{S}$